

Revised Details of SS Items

Pre-Bid meeting dated 17.06.16 to 18.06.2016

Tender Reference No: NPI-120310-EQP-S1-TD-24 Rev.00

Schedule No : 2


Schedule No – 2		
SI No	Item	Qty
1	Plain work table	153
	Vial inspection table	4
2	SS Ladders (3 Step)	15
3	SS Chairs	141
4	SS Bins	387
	SS Pedestal discard bin	
	SS Feed bins with Scoops (Lid type)	
	SS Garment discard bin	
5	SS Janitor closet	19
6	SS Lyophilizer loading table with bridge plate	1
7	SS Trolleys and Dress cabinets	116
	SS Platform trolley	
	Double Decker trolleys	
	Three Decker trolley	
	SS Openable trolley with doors	
	SS Dress cabinet/ used garment collection cabinet	
8	SS Perforated bins filling syringes*	0

*Included in schedule 1 sl no 5 & 6 as surgical bins

1	Work Tables
<p>Work tables are 2 types</p> <ul style="list-style-type: none"> • Plain work table • Vial visual Inspection Table 	
Plain Work Table :	
<ul style="list-style-type: none"> • Work tables shall be used in clean room areas • Work tables will not be having any shelves means plain table • MOC of the work table should be SS 304, Finish 240 Grit • 'C' Type support to the legs provided for sturdiness with a minimum four number of stiffener at top and bottom. • Up to the length of the 1500 mm 4 legs shall be considered and more than 1500 mm length table 6 legs shall be considered. • Vendor to consider proper support and stiffeners. • Top sheet of the table shall be of 16 Gauge thickness • 25 mm square pipe with 18 gauge thickness shall be consider for the fabrication. • For the legs Non- shedding Nylon bushes shall be considered. • Dimension as per the BOQ 	

Box

SS Work Tables				
Block	Item	Dimension	No's	Total
QC Block	SS Work Table	900 X 600 X 900 (LXWXH)	7	11
QC Block	SS Work Table	900 X 750 X 900 (LXWXH)	4	
BVF Block	SS Work Table	1200 X 600 X 900 (LXWXH)	6	27
QC Block	SS Work Table		1	
Rabies Block	SS Work Table		8	
SP Block	SS Work Table		2	
VVF Block	SS Work Table		6	
WH	SS Work Table		4	
VVF Block	SS Work Table	1200 X 750 X 900 (LXWXH)	1	66
Animal House Block	SS Work table		2	
SP Block	SS Work Table		1	
BCG Block	SS Work Table		25	
MBB Block (Hepatitis B Vaccine)	SS Work Table		12	
MBB Block (Hib Vaccine)	SS Work Table		23	
BVF Block	SS Work Table	1500 X 600 X 900 (LXWXH)	2	21
BVF Block	SS Work Table		6	
QC Block	SS Work Table		1	
Rabies Block	SS Work Table		3	
SP Block	SS Work Table		1	
VVF Block	SS Work Table		9	
WH	SS Work Table	1500 X 750 X 900 (LXWXH)	1	27
BCG Block	SS Work Table		3	
MBB Block (Hepatitis B Vaccine)	SS Work Table		11	
MBB Block (Hib Vaccine)	SS Work Table		1	
SP Block	SS Work Table		10	
VVF Block	SS Work Table		2	
MBB Block (Hepatitis B Vaccine)	SS Work Table	1800 X 750 X 900 (LXWXH)	1	1
Total Work Tables				153

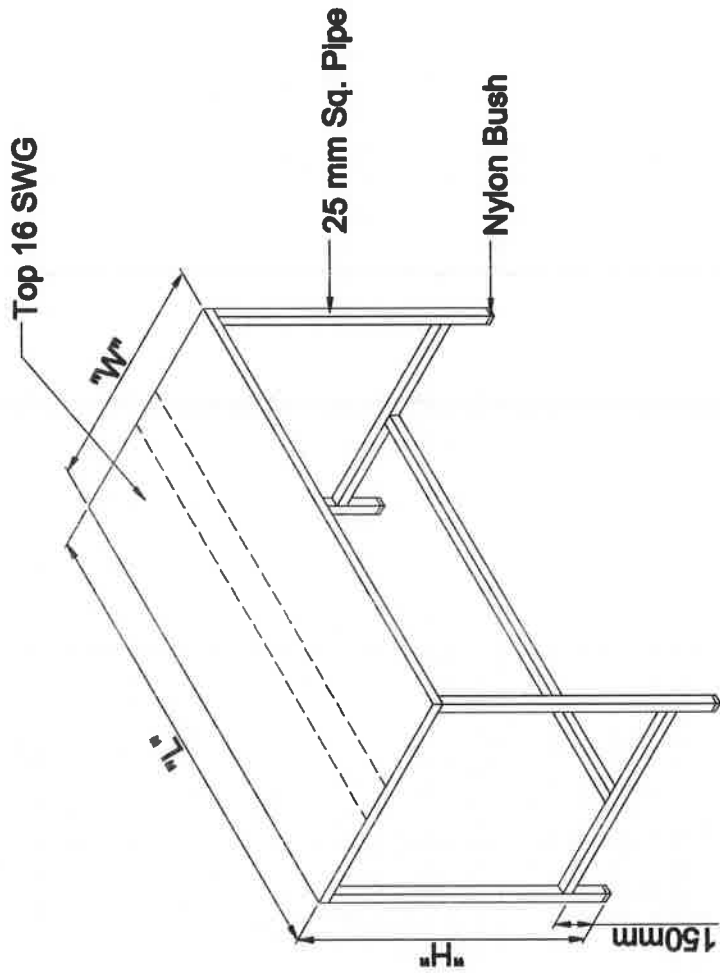

Item : SS Work Table
M.O.C.: SS 304 Qlty.* 16/18 Swg
Finish : Matt Finish


Drg. No. : 042

Vial Visual Inspection Table :

- Work tables shall be used in clean room areas for visual inspection of finished products.
- MOC of the work table should be SS 304, Finish 240 Grit
- 'C' Type support to the legs provided for sturdiness and total length of the 1220 mm with 4 legs shall be considered.
- Vendor should consider middle wall of SS Sheet with Black and White SS Vinyl sheet.
- Top sheet of the table shall be of 16 Gauge thickness and 40 mm square pipe thickness shall be consider for the fabrication.
- Top of platform shall be contained both side 6 nos rejection openings with 30mm Dia. The opening shall be covered with spring loaded cover.
- Total height of the table including SS canopy shall be 1450 (± 50 mm adjustable). The work table height shall be considered 850 mm(± 50 mm adjustable) from bottom.
- Both side 1 no of drawer with 6 compartment shall be consider for collection of rejected vials through rejection openings.
- Vendor should consider each side one no of 36 watt CFL tube light for visual inspection and ON/OFF switch also be considered as per diagram.
- For the legs Non- shedding Nylon bushes shall be considered.
- Vendor to consider proper support and stiffeners.
- Dimension as per the BOQ


BOQ

SS Work Tables - Vial Inspection				
Block	Item	Dimension	No's	Total
VVF Block	SS Work Table For Vial Inspection	1220 X 1050 - 1100 X 850 - 900 (LXWXH)	4	4
Total Vial Inspection tables				4

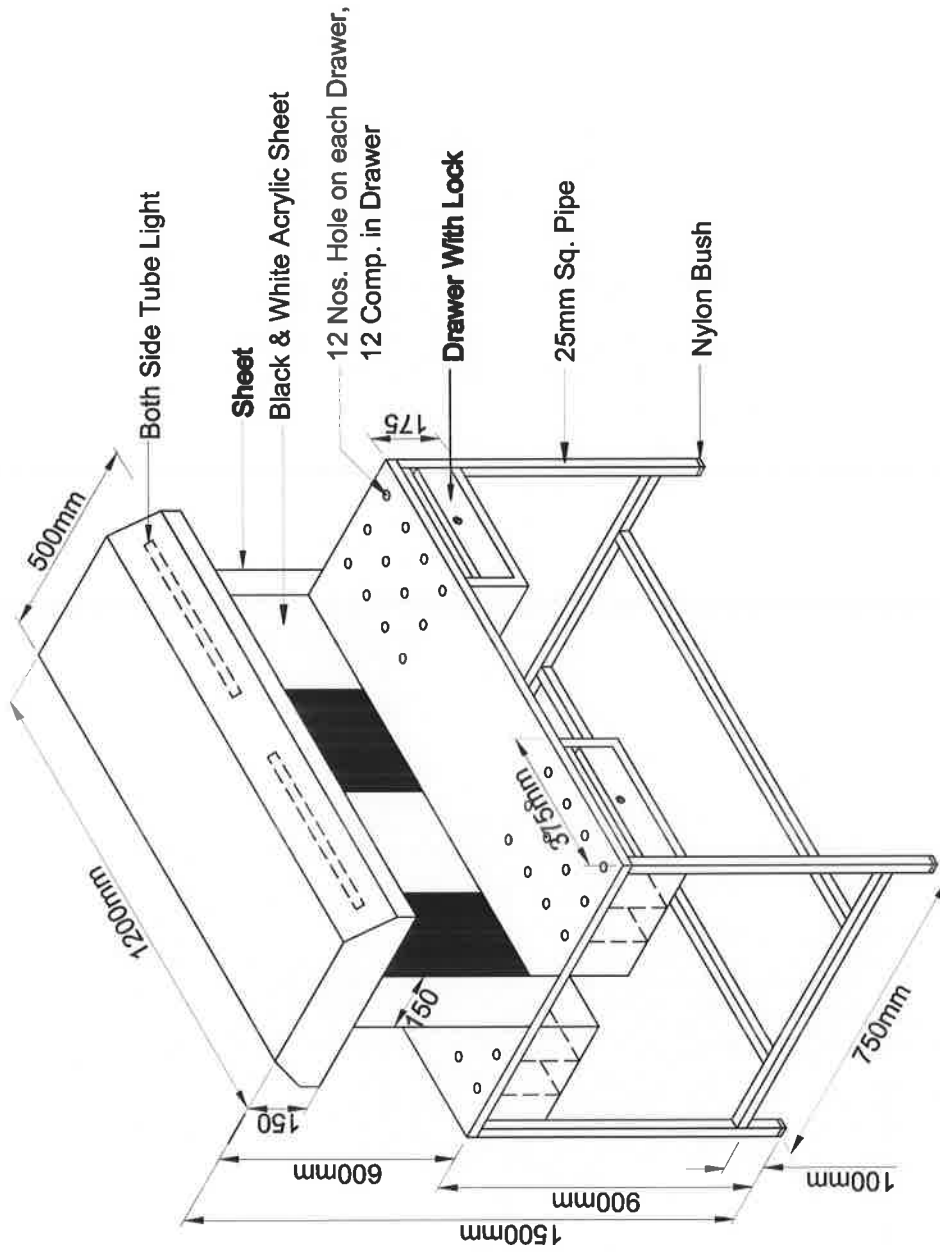

Item : SS Visual Inspection Table
M.O.C. : SS 304 Qlty.* 16/18 Swg
Finish : Matt Finish

VIAL VISUAL INSPECTION TABLE

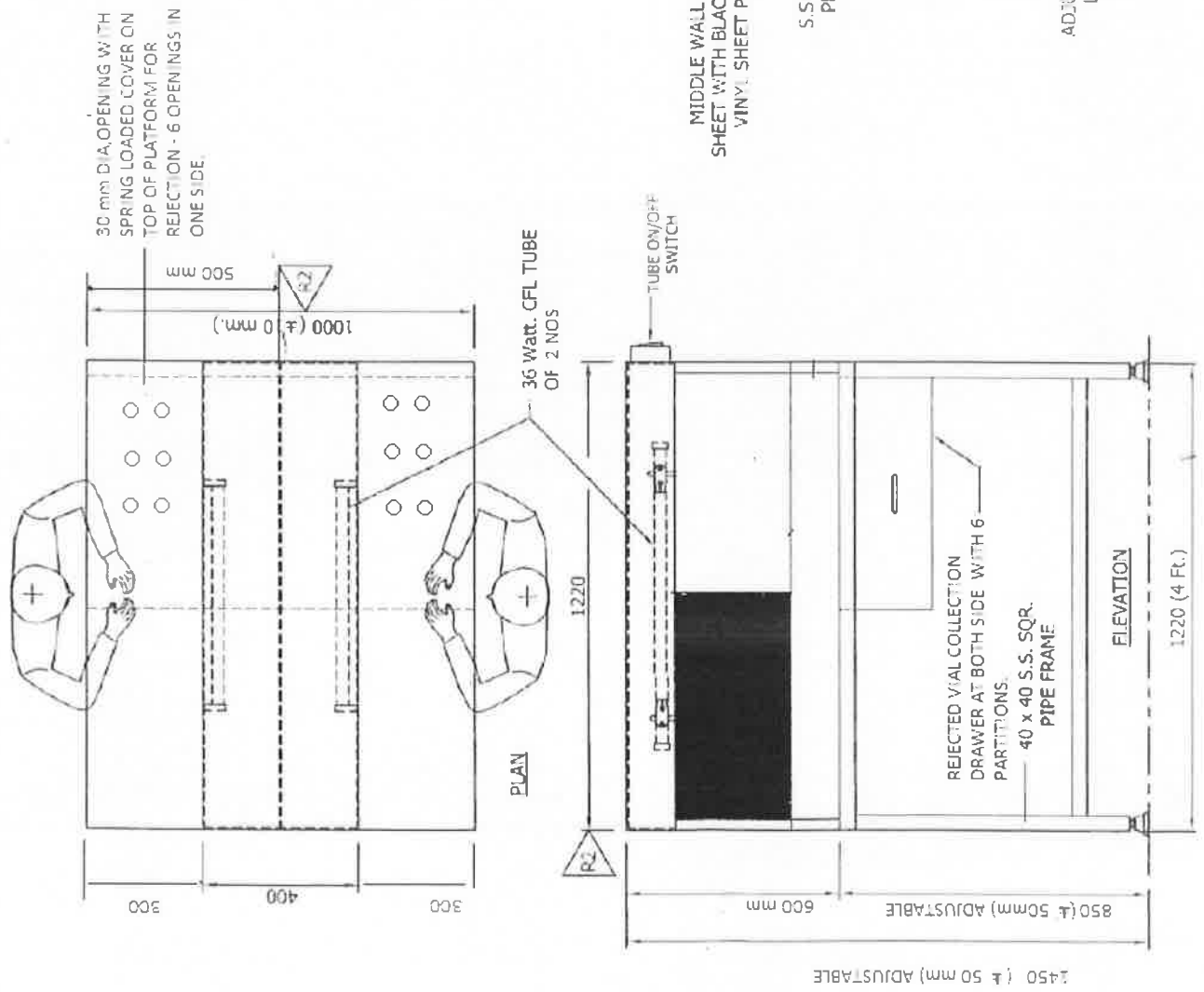


2	SS Ladders (3 Step) :
----------	------------------------------

- 3 Step ladders shall be used in clean rooms.
- MOC of the SS Ladders should be SS 304, Finish 240 Grit
- One Inch square pipe with 18 Gauge thickness should be used for fabrication.
- Total height of the Ladder shall be 900 mm, Length of SS Ladders 3 Steps shall be 600 mm Top Step (3rd) step of the ladder width shall be 300 mm and below two steps width shall be 250 mm.
- All sheets of the ladder steps shall be of dimple sheet with 16 Gauge thickness.
- Stiffeners shall be considered for sturdiness for the steps of the ladder.
- Ladders shall be matt finish and for the legs Non- shedding Nylon bushes shall be considered.

BOQ

SS Items - 3 Step Ladders				
Block	Item	Dimension	No's	Total
Animal House Block	SS 3 Step Ladders	800 X 600 X 900 (LXWXH)	3	15
BVF Block			1	
BCG Block			2	
Central Change Room and Laundry Block			1	
MBB Block (Hepatitis B Vaccine)			1	
MBB Block (Hib Vaccine)			1	
QC Block			1	
Rabies Block			2	
SP Block			1	
VVF Block			1	
WH			1	
Total SS 3 Step Ladders				15

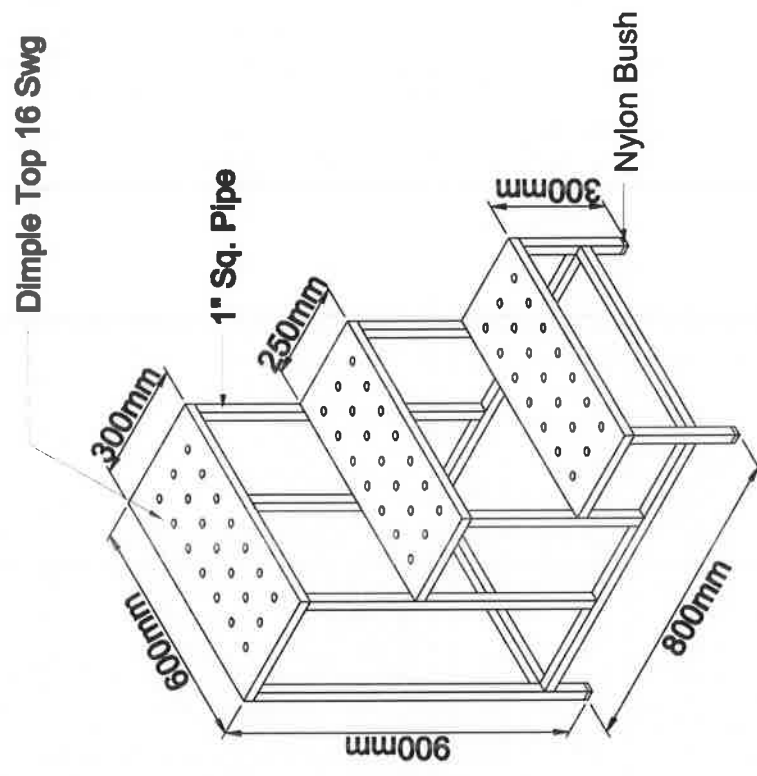

Item : SS Ladder 3 Steps
M.O.C. : SS 304 Qlty.* 16/18 Swg
Finish : Matt Finish


Drg. No. : 043

3 SS Chairs :

- SS Chairs shall be used in the Clean rooms.
- MOC of SS Chairs shall be of 304, Finish : 240 Grit
- 20 mm Dia pipe with 18 gauge thickness shall be used for fabrication.
- Back rest sheet and seating sheet shall be of 16 Gauge thickness .
- SS dia pipe shall be of mirror finish and back rest sheet and sitting sheet shall be of matt finish.
- Total height of the chair shall be 950 mm, Back rest sheet length shall be 450 mm with 200 mm Width, and seating sheet length shall be 450 mm with 320 mm width.
- Total seating height shall be 500 mm, from seating place to top of the back rest height shall be 450 mm, So that total height of the chair will be 950 mm.
- Stiffeners shall be considered for sturdiness
- For all the four legs interconnecting pipes shall be considered for sturdiness.
- Non- shedding Nylon bushes shall be considered for the legs.

BOQ

SS Chairs				
Block	Item	Dimension	No's	Total
Animal House Block	SS Chair	450 X 350 X 950 (LXWXH)	15	141
	SS Chair		4	
	SS Chair		4	
BCG Block	SS Chair		16	
	SS Chair		10	
MBB Block (Hepatitis B Vaccine)	SS Chair		25	
	SS Chair		26	
QC Block	SS Chair		19	
	SS Chair		14	
Rabies Block	SS Chair		8	
	SS Chair			
VVF RABIES				
VVF Measles				
WH				
Total Chairs				141

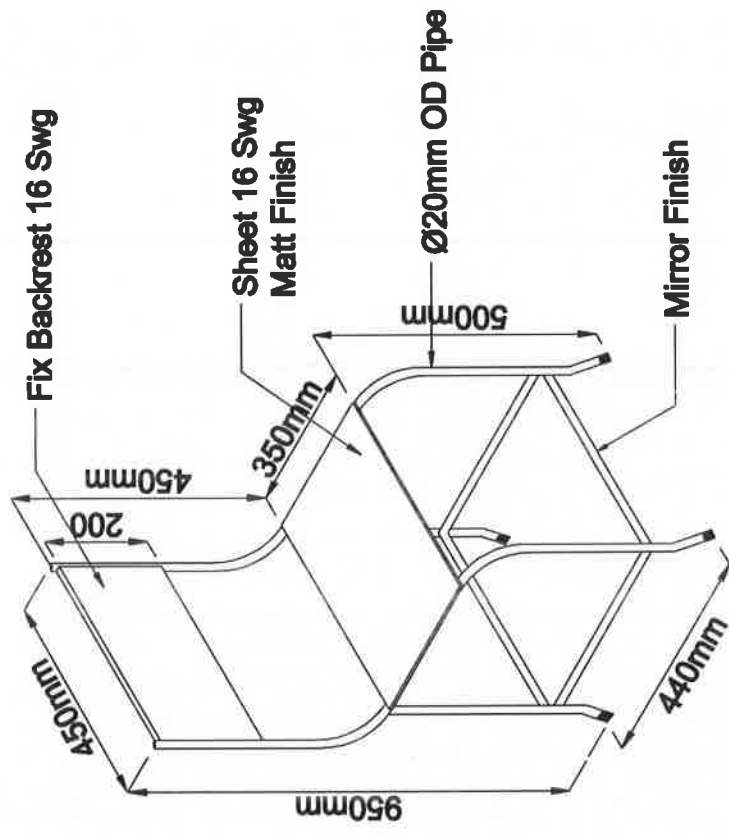

Item : SS Fix Chair

M.O.C. : SS 304 Qlty.* 16/18 Swg

Finish : Matt & Mirror Finish

Drg. No. 44


4	SS Bins :
<p>SS Bins are of different types</p> <ol style="list-style-type: none"> 1. SS Discard Bins 2. SS Feed bins with scoops (Lid Type) 3. SS Garment Discard Bin 	
SS Discard Bin :	
<ul style="list-style-type: none"> • Discard bins shall be of Pedestal type <p>Pedestal Discard Bin:</p> <ul style="list-style-type: none"> • MOC of Bins shall be SS 304, Finish : 240 Grit • Lid shall be opened by pressing with foot and handle shall be provided. • Smooth ground surface, with base hoop (Bottom sheet shall be beaded). • All the materials in foot operated dust bin shall be SS 304 • Sheet thickness shall be of 18 Gauge with Matt finish. • Dimension of the Foot operated Dust bin shall be 300 mm (Dia) X 450 mm (Height). 	
SS Feed bins with scoops (Lid Type) :	
<ul style="list-style-type: none"> • MOC of Bins shall be SS 304, Finish : 240 Grit • Lid shall be with top handle to remove or for placement on the bin. • Smooth ground surface, with base hoop (Bottom sheet shall be beaded) and 2 handles shall be provided. • All the materials Like Bin, Lid and Scoops shall be SS 304 • Sheet thickness shall be of 18 Gauge with Matt finish. • Scoop holding capacity shall be 0.5 Kg size • Dimension of the SS Feed bins shall be 300 mm (Dia) X 600 mm (Height). 	

SS Garment Discard Bin :

- MOC of SS Garment discard Bins shall be SS 304, Finish : 240 Grit
- SS Garment discard bin shall be in round shape.
- Top of the garment discard bin, Middle of the bin and bottom of the bin (Total 3 rims) shall be considered for the sturdiness.
- Top rim shall have 6 hooks from inside the bin this is to hang the lint free cloth or polythene cover for disposing the cloths.
- Four Legs with round shape pipes shall be considered for the support to the rims.
- 20 mm 18 gauge thickened pipes shall be considered for the Garment discard bin.
- To the base support of the legs non shedding nylon bushes shall be considered.
- Total dimension of the SS Garment discard shall be 450 (Dia) X 750 (Height)

BOQ

SS Bins					
Block	Item	Dimension	No's	Total	
Animal House Block	SS Pedestal Discard Bin	300 X 450 (Dia X H)	23	282	
	SS Pedestal Discard Bin		17		
	SS Pedestal Discard Bin		37		
Central Change Room and Laundry Block	SS Pedestal Discard Bin		18		
	SS Pedestal Discard Bin		27		
	SS Pedestal Discard Bin		23		
MBB Block (Hib Vaccine)	SS Pedestal Discard Bin		26		
	SS Pedestal Discard Bin		45		
	SS Pedestal Discard Bin		8		
VVF Block	SS Pedestal Discard Bin		49		
	WH		SS Pedestal Discard Bin		9
Total Pedestal Discard Bins				282	

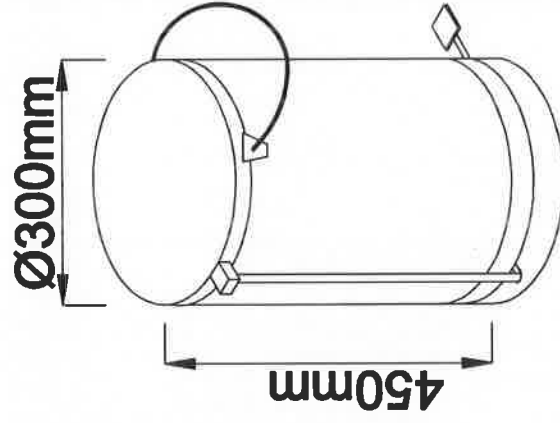


Item : SS Paddle Dustbin
M.O.C. : SS 304 Qlty. * 18 Swg
Finish : Matt Finish


Drg. No. : 69

B08

SS Bins				
Block	Item	Dimension	No's	Total
Animal House Block	SS Feed bins with Scoops (Lid Type)	300 X 600 (Dia x H)	12	12

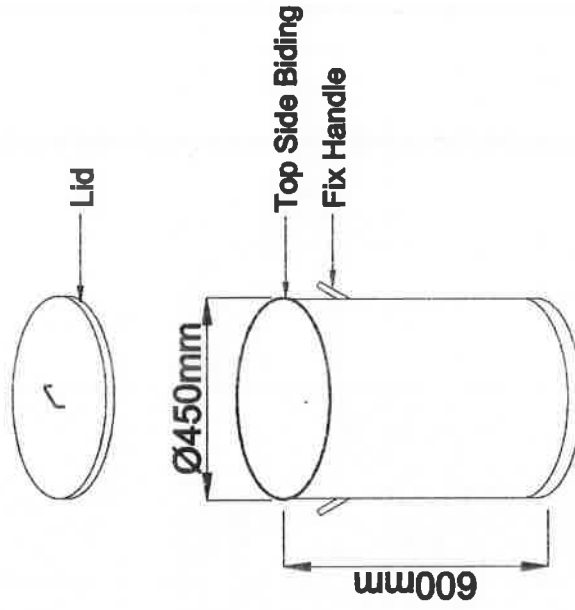

OPEN TYPE
500 Gms. Cap.

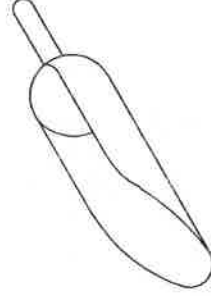


Item : SS Feed Bin With Scoop
M.O.C. : SS 304 Qty. * 18 Swg
Finish : Inside Mirror Finish /
Outside Matt Finish

Drg. No. 027

BOB

SS Bins				
Block	Item	Dimension	No's	Total
Animal House Block	SS Garment Discard Bin	450 X 750 (Dia X H)	10	93
	SS Garment Discard Bin		12	
	SS Garment Discard Bin		8	
BVF Block	SS Garment Discard Bin		6	
	SS Garment Discard Bin		7	
	SS Garment Discard Bin		7	
BCG Block	SS Garment Discard Bin		9	
	SS Garment Discard Bin		12	
	SS Garment Discard Bin		2	
Central Change Room and Laundry Block	SS Garment Discard Bin		16	
	SS Garment Discard Bin		4	
	SS Garment Discard Bin			
MBB Block (Hepatitis B Vaccine)				
MBB Block (Hib Vaccine)				
QC Block				
Rabies Block				
SP Block				
VVF Block				
WH				
Total Garment Discard Bins			93	

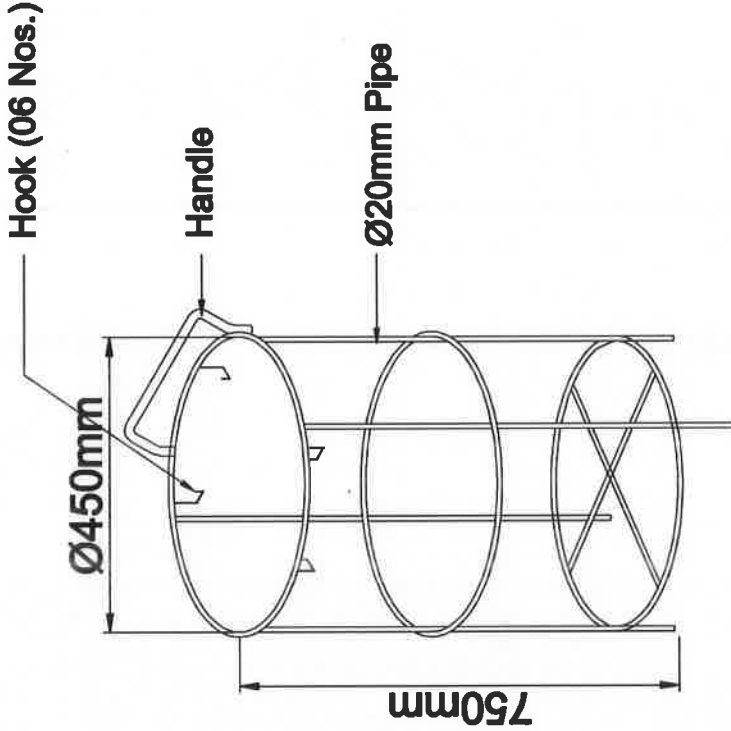

Item : SS Garment Discard Bin
M.O.C. : SS 304 Qlty. * 18 Swg
Finish : Matt Finish


Drg. No. : 70

5 SS Janitor closet :

- SS Janitor shall be used to keep the cleaning materials like brooms, mops, mop sticks and cleaning agents or disinfectants in clean rooms of dedicated janitor rooms.
- In the Janitor at the centre one partition shall be made, Left side partition shall be plain (means without shelves) and right side partition shall be with 3 shelves each at the height of 450 mm so total compartments will be 4.
- Two partitions of the Janitor shall have two doors with lock and key.
- In the left side of the partition on top hooks shall be considered for hanging mops and mop sticks.
- MOC of the Janitor closet shall be SS 304, Finish 240 Grit .
- SS Sheet shall be of 18 Gauge thickened, SS 304 with 240 Grit finish.
- Top portion of the Janitor closet shall be in a slope position to avoid any deposition and for easy cleaning.
- Dimensions of the janitor closet shall be 900mm (Length) X 450mm (Width) X 1800 mm (Height).
- Due to slope on the top portion Janitor back side height shall be 1900 mm and at the front it shall be 1800 mm.
- Legs shall be provided with non-shedding nylon bushes at the bottom with 150 mm clearance from the base.
- Vendor to consider proper support and stiffeners.

BoQ

SS Janitor closet				
Block	Item	Dimension	No's	Total
Rabies Block	Janitor Closet	900 X 450 X 1800 (LXWXH)	4	4
Animal House Block	Janitor Closet	1000 X 450 X 1800 (LXWXH)	1	1
SP Block	Janitor Closet	1200 X 450 X 1800 (LXWXH)	1	
VVF Block	Janitor Closet	1200 X 450 X 1800 (LXWXH)	1	
VVF Block	Janitor Closet	1200 X 450 X 1800 (LXWXH)	2	
Animal House Block	Janitor Closet	1200 X 450 X 1800 (LXWXH)	4	
BVF Block	Janitor Closet	1200 X 450 X 1800 (LXWXH)	2	
BCG Block	Janitor Closet	1200 X 450 X 1800 (LXWXH)	1	
MBB Block (Hepatitis B Vaccine)	Janitor Closet	1200 X 450 X 1800 (LXWXH)	1	
MBB Block (Hib Vaccine)	Janitor Closet	1200 X 450 X 1800 (LXWXH)	1	
QC Block	Janitor Closet	1200 X 450 X 1800 (LXWXH)	1	14
			Total Janitor closet	19

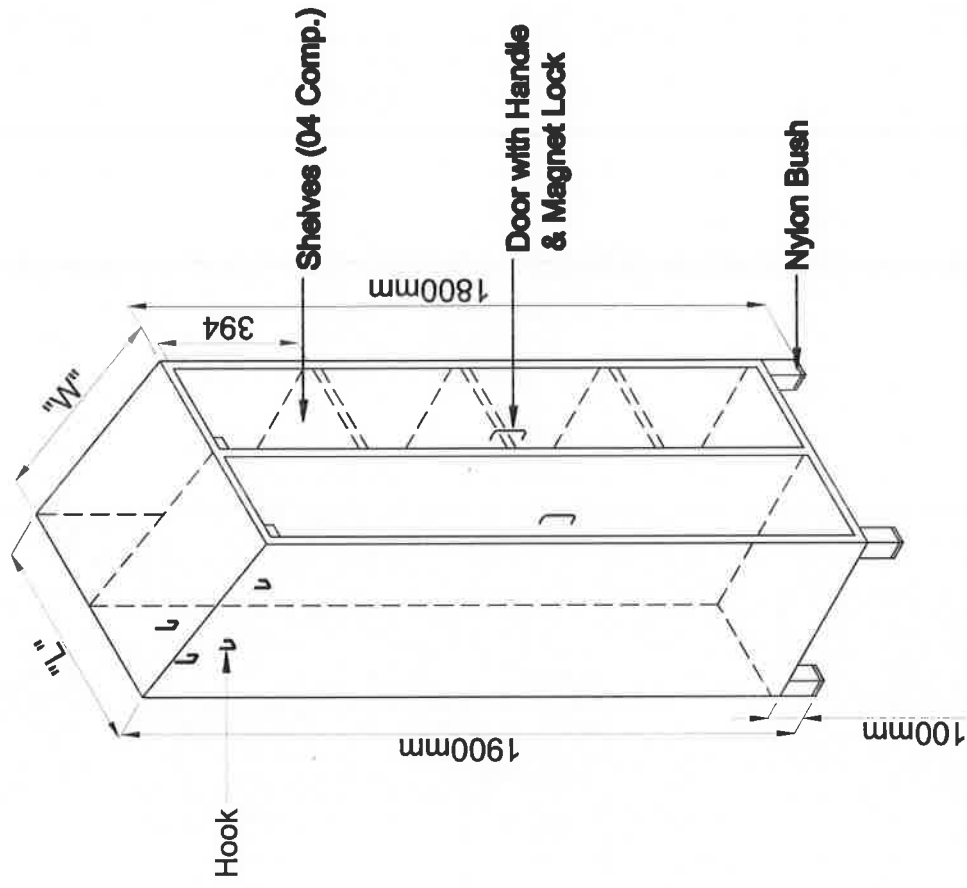

Item : SS Housed Keeping Cupboard
M.O.C. : SS 304 Qlty.* 18 Swg
Finish : Matt Finish

Drg. No. 45


6 SS Lyophilizer loading table with bridge plate :

- Lyophilizer loading table shall be used for keeping the Lyophilizer on the top of the Table.
- Lyophilizer loading table shall be of SS Height adjustable type.
- Folding type tray shall be provided with stiffener , dimension shall be of 904 mm (L) X 800 mm (Width)
- Total size of the table shall be 2000 mm (L) X 990 mm (W) X 900 mm to 1000 mm (Height adjustable).
- 38 mm square pipe with 18 Gauge thickness should be used for the fabrication.
- Top sheet of the table shall be with 16 Gauge thickness.
- MOC of the Lyophilizer loading table should be SS 304, Finish 240 Grit
- Vendor to consider proper support and stiffeners.
- Table shall be of Matt finish.
- Dimension as per the BOQ

BoQ

SS Items - Lyo loading table with bridge plate				
Block	Item	Dimension	No's	Total
BCG	Lyo loading table with bridge plate	2000 X 990 X 900 - 1000 mm (LXWXH) Height adjustable	1	1
				Folding plate size shall be 904 mm X 800 mm)

Hydroloading Table with Birdge Plate

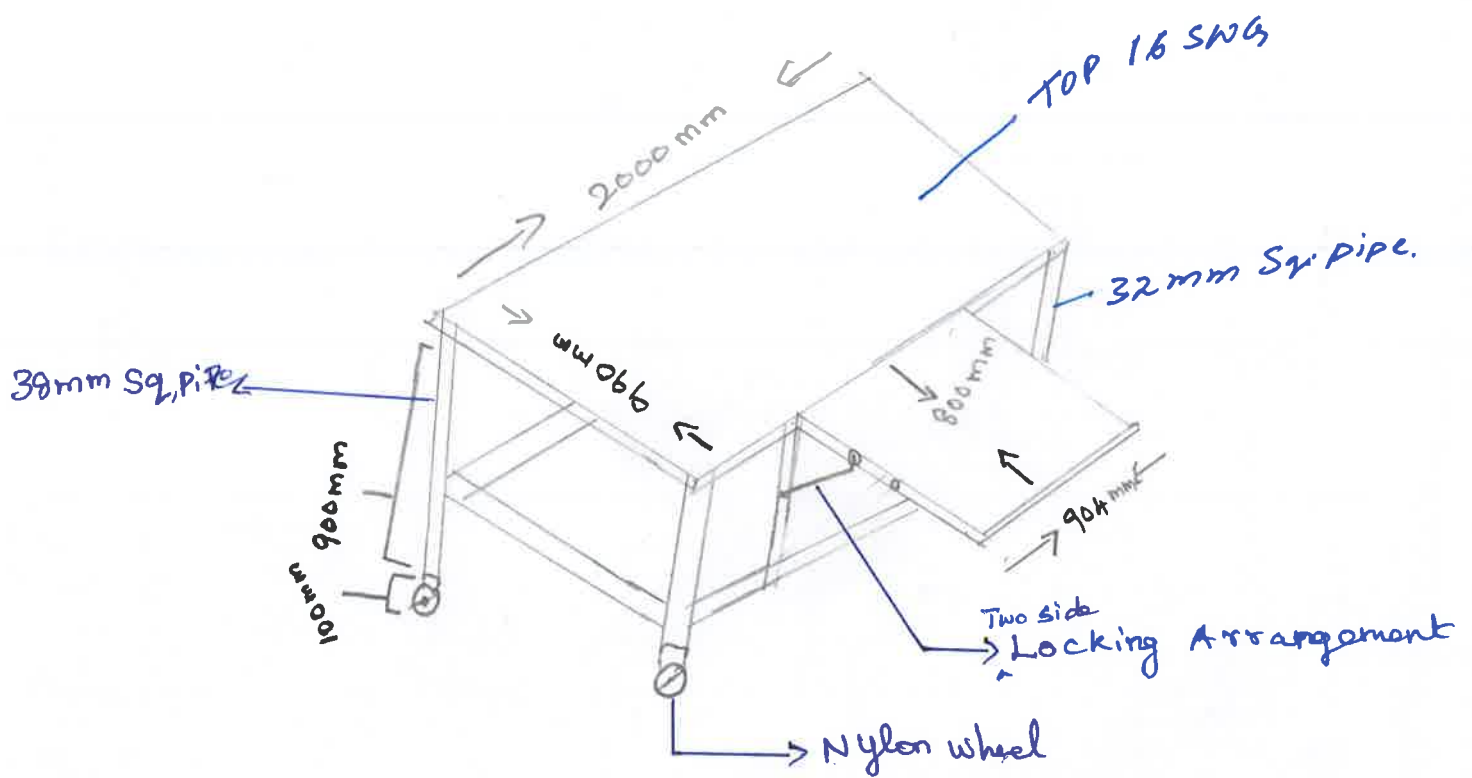

900 to 1000 mm adjustable height

7	SS Trolleys and Dress cabinets :
---	---

SS Trolleys are of different types

- Platform trolleys
- Double Decker trolleys
- Three Decker trolleys
- Openable (Closed with doors) trolley

SS Trolleys shall be used in the Clean rooms for transfer of the materials from one area to other area.

Platform Trolleys :

- Some of the trolleys shall be with platform edges (4 side of the platform trolley with guard shall be considered up to the height of 100 mm) this is to prevent bottle falling down from the shelves. Some of the platform trolleys shall be of Plain type. (Plain surface).
- No of guarded platform trolleys and plain platform trolleys with sizes are mentioned in the BOQ.
- MOC of the trolleys shall be SS 304, Finish : 240 Grit
- Platform or base sheet of the trolleys should be 16 Gauge thickness with stiffener for strength & load bearing capacity.
- Handle should be made of 1 Inch dia pipe with 16 Gauge thickness with stiffener in between the rod for firm grip & load pushing capacity.
- SS Platform with heavy duty PU wheels of size 4 Inch (Dia) X 2 Inch (Width) should be considered.
- Two wheels shall be with free movement and other two wheels shall be with brakes.
- Dimension of the platform trolleys shall be 800 mm (L) X 650 mm (W) X 900 mm (Handle Height).
- Platform trolley shall be of Matt finish and clamps of the wheels shall be of SS 304.

Double Decker Trolleys :

- Double Decker trolleys shall be with edges (4 side of the two decker or tray guard shall be considered up to the height of 50 mm) this is to prevent bottle falling down from the shelves.
- No of Double Decker trolleys with sizes are mentioned in the BOQ.
- MOC of the trolleys shall be SS 304, Finish : 240 Grit
- SS Sheet for the trolleys should be 16 Gauge thickness with stiffener for strength & load bearing capacity.
- Handle should be made of 1 Inch dia pipe with 16 Gauge thickness with stiffener in between the rod for firm grip & load pushing capacity.
- SS Platform with heavy duty PU wheels of size 4 Inch (Dia) X 2 Inch (Width) should be considered.
- Two wheels shall be with free movement (Revolving type) and other two wheels shall be with brakes.
- Dimension of the Double Decker rolleys shall be 800 mm (L) X 650 mm (W) X 900 mm (Handle Height).
- Trolley shall be of Matt finish and clamps of the wheels shall be of SS 304.

Triple Decker (or) 3 Tray Trolleys :

- Triple Decker trolleys (or) 3 Tray trolley shall be with edges (4 sides of the three decker or tray guard shall be considered up to the height of 50 mm) this is to prevent bottle falling down from the shelves.
- No of Three Decker (or) 3 Tray trolleys with sizes are mentioned in the BOQ.
- In between the tray height shall be 300 mm.
- MOC of the trolleys shall be SS 304, Finish : 240 Grit
- SS Sheet for the trolleys should be 16 Gauge thickness with stiffener for strength & load bearing capacity.
- Handle should be made of 1 Inch dia pipe with 16 Gauge thickness with stiffener in between the rod for firm grip & load pushing capacity.
- SS Platform with heavy duty PU wheels of size 4 Inch (Dia) X 2 Inch (Width) should be considered.
- Two wheels shall be with free movement (Revolving type) and other two wheels shall be with brakes.
- Dimension of the Double Decker rolleys shall be 800 mm (L) X 650 mm (W) X 1100 mm (Handle Height).
- Trolley shall be of Matt finish and clamps of the wheels shall be of SS 304.

Openable (Closed with doors) trolley :


- Openable Trolley shall be with doors to carry the carrier boxes from one area to another area.
- 2 Shelves shall be provided at equally distributed height in the closed type trolley so that total compartments could be three. Total trolley shall be closed from all sides.
- No of Openable trolleys with sizes are mentioned in the BOQ.
- MOC of the trolleys shall be SS 304, Finish : 240 Grit
- Two doors with handles and drop lockers shall be considered at the front.
- Partition is not required in the compartments. So compartments shall be of Plain.
- SS Sheet for the Openable trolley should be 18 Gauge thickness with stiffener for strength & load bearing capacity.
- Handle shall be made of 1 Inch dia pipe with 18 Gauge thickness for firm grip & load pushing capacity.
- SS Platform with heavy duty PU wheels of size 4 Inch (Dia) X 2 Inch (Width) should be considered.
- Two wheels shall be with free movement (Revolving type) and other two wheels shall be with brakes.
- Dimension of the Openable Trolley shall be 900 mm (L) X 900 mm (W) X 1300 mm and Handle height shall be at 900 mm on both sides.
- Trolley shall be of Matt finish and clamps of the wheels shall be of SS 304.

SS Dress Cabinets / Used garments collecting cabinets :

- SS Dress cabinets shall be placed in the central change room block for to transfer the used garments from Different block to the Laundry area.
- Front side of the Dress cabinet shall be closed completely and on top side of the front sheet cut out shall be provided to drop the used clothes. To the cut out swing door shall be provided from inside of the cabinet. Cabinet cut out and door size shall be 600 (L) X 400 (W).
- Back side of the cabinet door shall be provided to collect the used garments for washing purpose. Back side door shall have the view panel with toughened glass approx. dimension of view panel shall be 450 (L) X 450 (W) at centre of the door.
- MOC of the Dress cabinets shall be SS 304, Finish : 240 Grit
- Back side door should be provided with magnetic latches and handle with lock and key.
- All sides of the dress cabinet shall be of 18 Gauge thickness.
- Bottom sheet and back side door shall be of 16 Gauge sheet.
- Stiffeners shall be considered for the sturdiness and strengthening of the Dress cabinets.
- Dimension of the SS Dress cabinets shall be 1000 mm (L) X 1000 mm (W) X 900 mm .
- SS Dress cabinets shall be of Matt finish


BOQ

SS Plat form Trolleys							
Block	Item	Dimension	No's	Total	Specification	Remarks	
BVF Block	SS Trolley	800 X 650 (LXW) 900 Handle height	4	41	Platform Trolleys	With Guard	
	SS Trolley		4			With Guard	
	SS Trolley		6			With Guard	
	SS trolley		3			With Guard	
Rabies Block	SS Trolley		2			With Guard	
	SS Trolley		4			With Guard	
VVF Block	SS Trolley		5			Without Guard	
	SS Trolley		2			Without Guard	
Animal House Block	SS Trolley		1			Without Guard	
	SS Trolley		10			Without Guard	
Central Change Room and Laundry Block	SS Trolley						
WH	SS Trolley						
SP Block	SS Trolley						
Total Platform Trolleys				41			

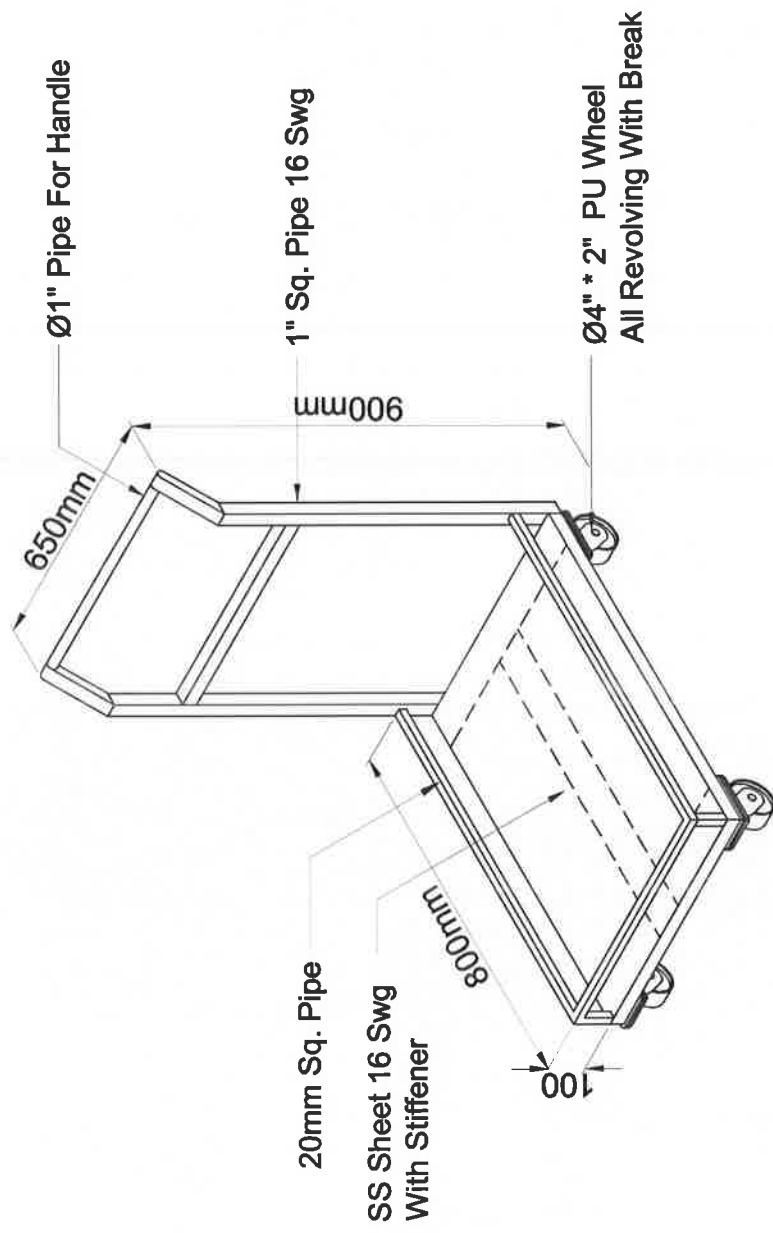

Drg. No. 66

Item : SS Platform Trolley
M.O.C. : SS 304 Qlty* 16 Swg
Finish : Matt Finish


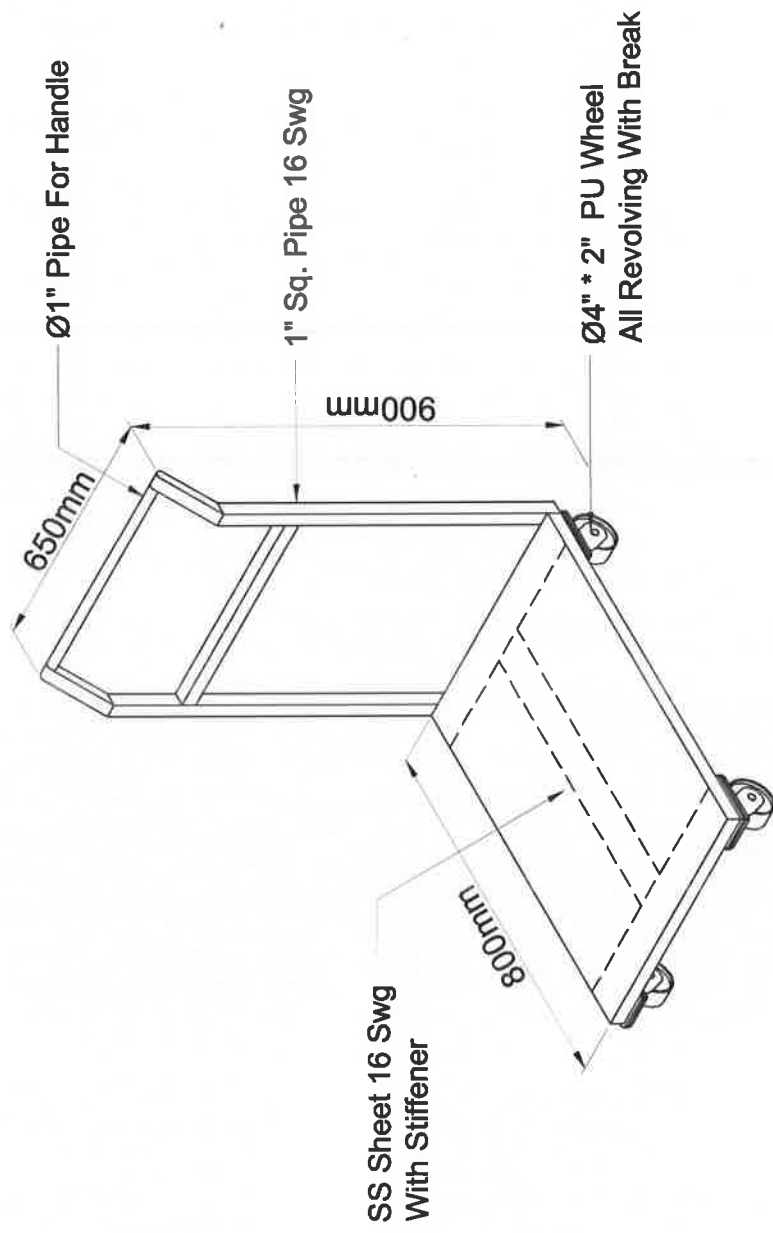


Item : SS Platform Trolley
M.O.C. : SS 304 Qlty* 16 Swg
Finish : Matt Finish

Drg. No. 67

B02



SS Double Decker Trolleys						
Block	Item	Dimension	No's	Total	Specification	Remarks
BVF Block	SS Trolley	800 X 650 (LXW) 900 Handle Height	4	56	Double Decker	With Guard
BCG Block	SS Trolley		14			With Guard
MBB Block (Hepatitis B Vaccine)	SS Trolley		3			With Guard
	SS trolley		11			With Guard
Rabies Block	SS Trolley		8			With Guard
VVF Block	SS Trolley		16			With Guard
Total Double Decker Trolleys				56		

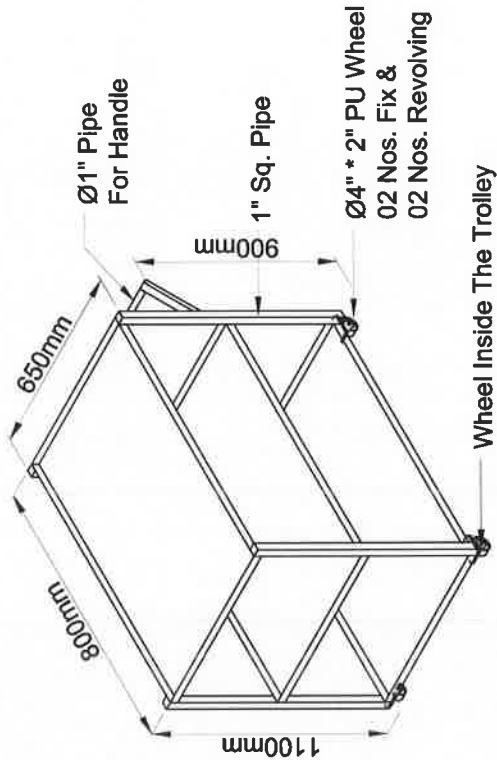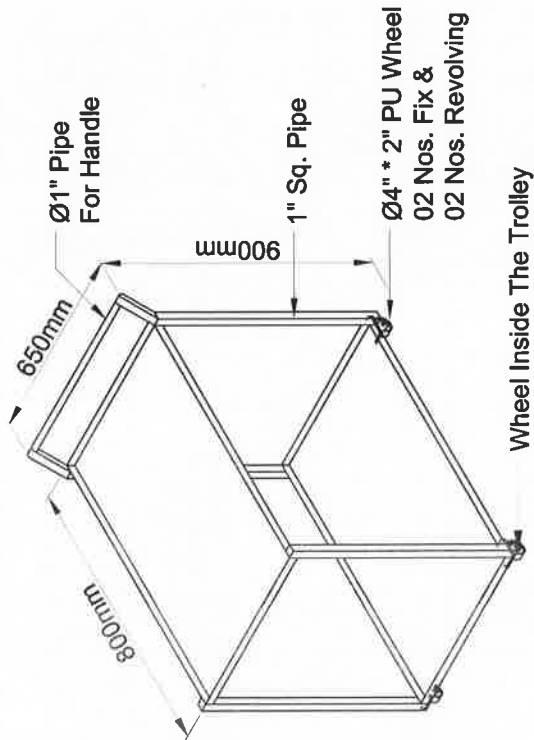

Item : SS Two Tier Trolley & Three Tier Trolley
M.O.C. : SS 304 Qlty. * 16 Swg
Finish : Matt Finish

BoQ


SS Three decker trolleys						
Block	Item	Dimension	No's	Total	Specification	Remarks
Animal House Block	SS Trolley	800 X 650 (LXW) 1100 Handle height	13	13	Three Decker	3 Tray trolley with borders on all 4 sides Tray to tray ht would be 450
Total Three decker trolleys				13		

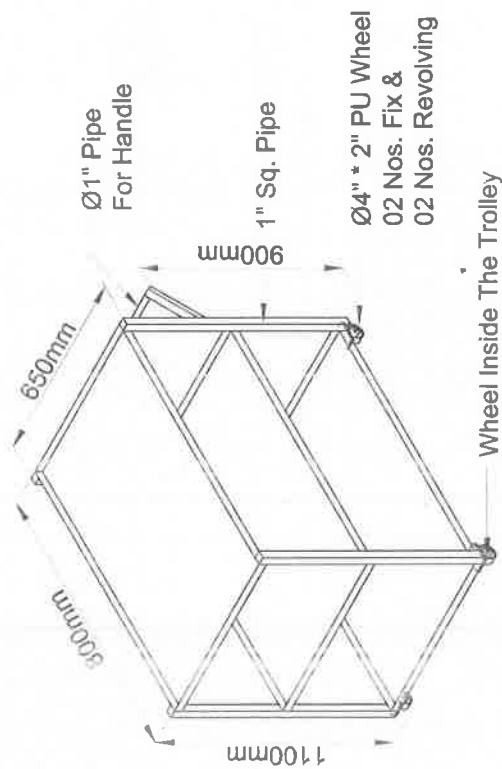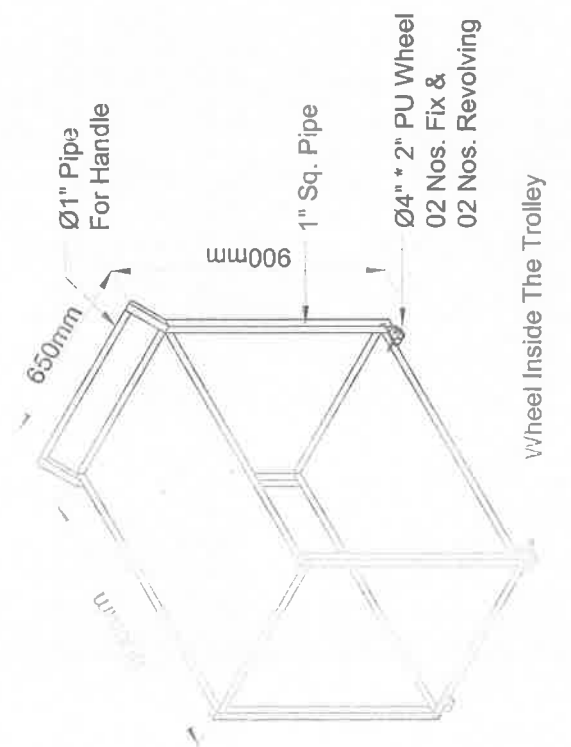

Item : SS Two Tier Trolley & Three Tier Trolley
M.O.C. : SS 304 Qlty. * 16 Swg
Finish : Matt Finish

BOQ

SS Openable Trolleys with doors					
Block	Item	Dimension	No's	Total	Specification
VVF Block	SS Trolley	900 X 900 X 1300 (LXWXH) 900 Handle Height	2	6	Openable Box type
MBB (Hib)	SS Trolley		1		
QC	SS Trolley		2		
BCG Block	SS Trolley		1		
Total Oenable trolleys with doors				6	
Completely closed Trolley with Doors					

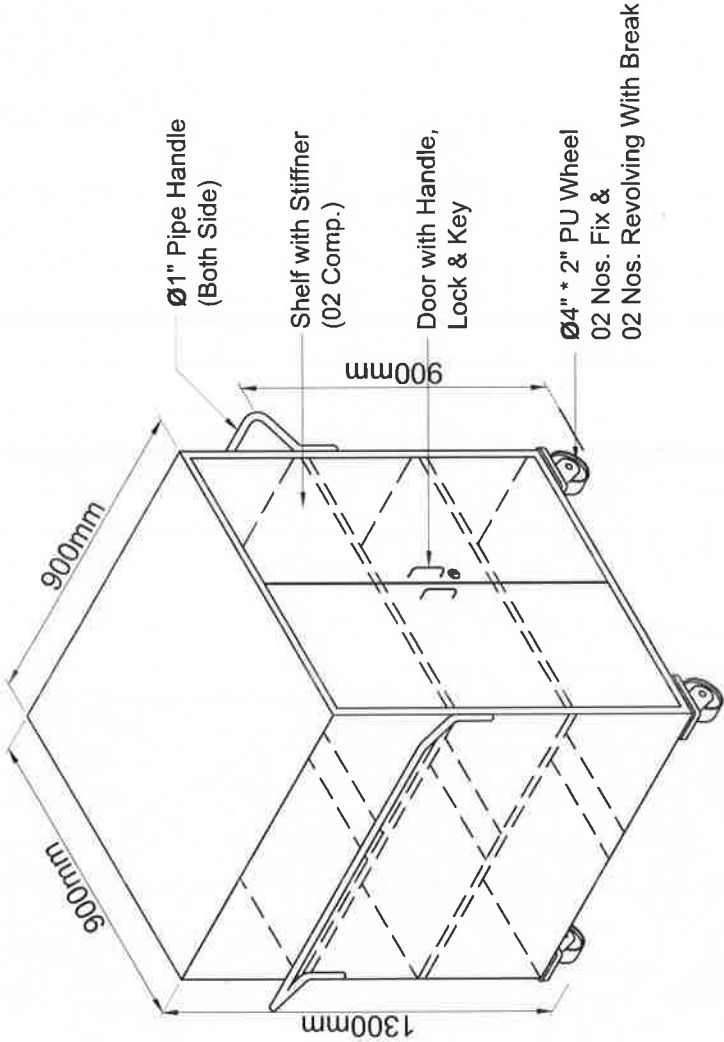


SS Surgical Bin / Box Storage Trolley

Item : SS Trolley (Openable Box Type)
M.O.C. : SS 304 Qlty* 18 Swg
Finish : Matt Finish

Drg. No. 68

Box

SS Dress Cabinets / Used garments collecting cabinets					
Block	Item	Dimension	No's	Total	Specification
Central Change Room and Laundry Block	SS Used Dress Collecting Cabinets	1000 X 1000 X900 (LXWXH)	10	10	Closed type with vie panel door on back side and front side flap door
Total SS Dress Cabinets / Used garments collecting cabinets				10	

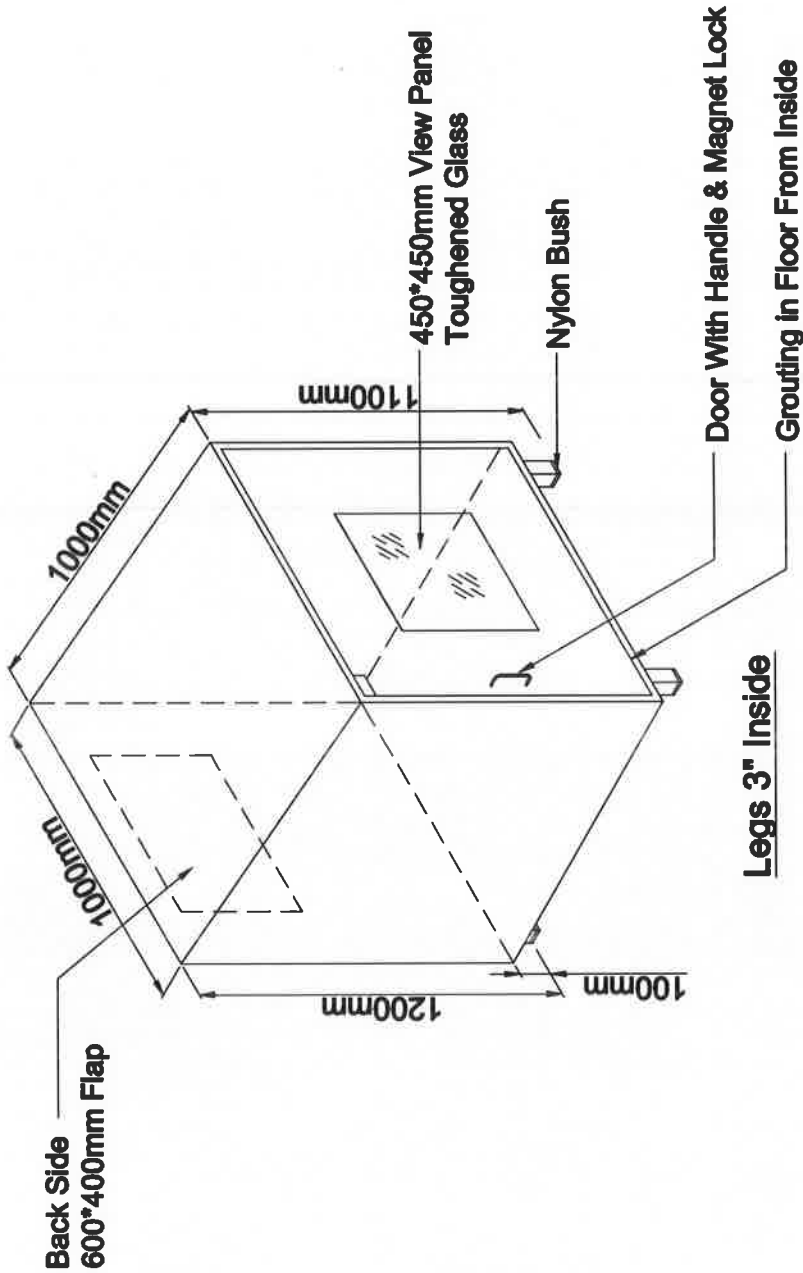

Item : SS Used Dress Collection Cabinet
M.O.C. : SS 304 Qty. * 16/18 Swg
Finish : Matt Finish

Drg. No. : 71