


HLL LIFECARE LIMITED, Chennai

Revival of BCG Vaccine Laboratory, Guindy, Chennai

HLL pharmaplan®	User Requirement Specifications				
	Equipment/System	Decontamination - Autoclave			
	Identification	BF-DAT 01 & FG-DAT 01	Document	URS/BF-DAT01_FG-DAT 01	
	Effective Date	2013.03.04	Revision	02	

User Requirement Specifications Decontamination - Autoclave


Equipment ID	BF-DAT 01
	FG-DAT 01

Revision index

Revision	Date	Reason for revision
00	12.06.2012	First Draft for Client's review
01	2013-01-03	Format Has been changed by HLL
02	2013.02.07	Client's comment (dated 2013.01.07) incorporated

HLL LIFECARE LIMITED, Chennai

Revival of BCG Vaccine Laboratory, Guindy, Chennai

nne pharmaplan®	User Requirement Specifications				
	Equipment/System	Decontamination - Autoclave			
	Identification	BF-DAT 01 & FG-DAT 01	Document	URS/BF-DAT01_FG-DAT 01	
	Effective Date	2013.03.04	Revision	02	

URS Annexure List

URS Annex No.	Detail
1.	Layouts showing location of the Decontamination Autoclave in bulk decontamination area and formulation washing area
2.	List of Preferred Make of components
3.	List of recommended spares

HLL LIFECARE LIMITED, Chennai

Revival of BCG Vaccine Laboratory, Guindy, Chennai



nne pharmaplan®	User Requirement Specifications				
	Equipment/System	Decontamination - Autoclave			
	Identification	BF-DAT 01 & FG-DAT 01	Document	URS/BF-DAT01_FG-DAT 01	
	Effective Date	2013.03.04	Revision	02	

Table of Contents

1.0	APPROVAL SIGNATURE	4
2.0	EQUIPMENT DESCRIPTION	5
3.0	PROCESS DESCRIPTION	7
3.1	INPUT & CHARGING METHOD	7
3.2	BRIEF PROCESS STEPS	7
3.3	OUTPUT & DISCHARGING METHOD	8
4.0	PRODUCTIVITY REQUIREMENT	8
4.1	DESIRED/ SUGGESTED CAPACITY	8
4.2	STANDARD BATCH SIZE	8
4.3	CHANGE OVER TIME	8
4.4	OTHER PRODUCTIVITY REQUIREMENT	8
5.0	CONTAINMENT	8
6.0	GMP REQUIREMENTS	8
6.1	PROCESS CONTROL	8
6.2	FAILURE MODE DETECTION	10
6.3	IN –PROCESS CONTROL	11
6.4	LEVEL OF INSTRUMENTATION	11
6.5	BATCH DATA DISPLAY AND RECORD PRINTING	12
6.6	GMP REQUIREMENTS (OTHERS)	12
6.7	SPECIFIC REQUIREMENTS	13
7.0	CONSTRAINTS	14
7.1	EQUIPMENT LOCATION AND AVAILABLE SPACE	14
7.2	UTILITY REQUIREMENT	16
8.0	ABBREVIATION	17

HLL LIFECARE LIMITED, Chennai

Revival of BCG Vaccine Laboratory, Guindy, Chennai

nne pharmaplan®	User Requirement Specifications				
	Equipment/System	Decontamination - Autoclave			
	Identification	BF-DAT 01 & FG-DAT 01	Document	URS/BF-DAT01_FG-DAT 01	
	Effective Date	2013.03.04	Revision	02	

1.0 Approval Signature

This document is prepared by the Validation and GMP compliance team of “NNE Pharmaplan India for the project “Revival of BCG Vaccine Laboratory” (**project number:-110729**) of BCG Vaccine Laboratory, Guindy, Chennai under the authority of their Project Manager. Hence, this document before being effective shall be approved by the QA team and authorized by the appropriate Project Authority.

Prepared by		
Name/ Designation	Signature	Date
Ms. Sandhya Samant Sr. Engineer – Projects (Biotech) NNE Pharmaplan India Ltd.		


Checked by		
Name/ Designation	Signature	Date
Dr. Naveen Nagaraj Sr. Project Manager - Head Biotech Division NNE Pharmaplan India Ltd.		
Mr. Arun Nagendra Manager-Validation and GMP NNE Pharmaplan India Ltd.		

Approved by		
Name/ Designation	Signature	Date
Mr. Narendra Prasad Director-Technical NNE Pharmaplan India Ltd		
HLL Lifecare Limited		

Authorized by		
Name/ Designation	Signature	Date
Project Authority BCG Vaccine Laboratory		

HLL LIFECARE LIMITED, Chennai

Revival of BCG Vaccine Laboratory, Guindy, Chennai

nne pharmaplan®	User Requirement Specifications				
	Equipment/System	Decontamination - Autoclave			
	Identification	BF-DAT 01 & FG-DAT 01	Document	URS/BF-DAT01_FG-DAT 01	
	Effective Date	2013.03.04	Revision	02	

2.0 Equipment description

Autoclave shall be used to decontaminate various items that come in contact with organisms and material according Biosafety class 2.


S.No	Identification no.	Purpose	Chamber inner Dimension (W x H x D in mm)
1.	BF-DAT 01	tubes, tuberlin flask, conical flask, flat round bottom flasks	600mm x 600mm x 600 mm
2.	FG-DAT 01	15L bulk bottle; buffer vessel (0.5L) along with nozzle arrangement and flexible tubing	600mm x 600mm x 600mm

- Operation programs for liquid, solid and porous goods and biohazard waste.
- Package Unit including the following:
 - Autoclave chamber
 - Supporting structure
 - Integrated vacuum system
 - filters for compressed air and for exhaust air
 - Piping (valves, safety devices, filters, steam traps, pipes, fittings, etc.)
 - Pressure reducing valve in plant steam line
 - All mating flanges/fittings, gaskets, bolts and screws for utility supplies, returns and drains
 - Instrumentation
 - Control System with printer (able to give trends in color)

Design, function and control of the units has to be **GMP compliant**

HLL LIFECARE LIMITED, Chennai


Revival of BCG Vaccine Laboratory, Guindy, Chennai

nne pharmaplan®	User Requirement Specifications				
	Equipment/System	Decontamination - Autoclave			
	Identification	BF-DAT 01 & FG-DAT 01	Document	URS/BF-DAT01_FG-DAT 01	
	Effective Date	2013.03.04	Revision	02	

Note:


I.	This Technical Specification is the basis for an inquiry to a vendor and therefore the basis for the vendor's proposal.
II.	The vendor is asked to state in "REMARKS" column with "yes" if the described requirement will be completely fulfilled and with "no" in case the requirement will not or cannot be fulfilled with the proposed equipment. In case of a deviation a comment must be inserted or enclosed as a separate annexure by referring to the respective URS specification number.
III.	The vendor must clearly comment each item of the Technical Specification. The comments must be in English language. If extra cost for necessary options becomes necessary the item must be clearly stated.
IV.	In case that the requirement includes a question or request or information from the vendor, the answer / information should be stated in the "REMARKS" column.
V.	The final version of this document including the vendor's comments will become basis of a potential purchase order or contract.
VI.	The Technical Specification serves to define a summary of all vendor's requirements concerning scope of delivery and services.
VII.	The vendor is responsible for technically unobjectionable function of the equipment. This TS is not intended to dictate a technical design to the vendor. If agreed upon with the vendor, the vendor can apply his practically proven design.
VIII.	<p>Special Instruction</p> <p>a. If no comments against any specification shall be considered as "NO" and</p> <p>b. If there is no reply / comments against the complete URS by the vendor then it shall be treated as unresponsive / technically non compliant and rejected.</p>
IX.	Refer document Installation Requirement Specification and Specific Instructions with URS; NPI_110729_IRS_BCG_01
X.	Refer Tender document with URS; NPI/110729/EQP/TD/05

HLL LIFECARE LIMITED, Chennai				
Revival of BCG Vaccine Laboratory, Guindy, Chennai				
HLL Pharmaplan®	User Requirement Specifications			
	Equipment/System	Decontamination - Autoclave		
	Identification	BF-DAT 01 & FG-DAT 01	Document	URS/BF-DAT01_FG-DAT 01
	Effective Date	2013.03.04	Revision	02
Specifications				Remarks
3.0 Process Description				
3.1 Input & Charging method				
3.1.1 The autoclave shall be suitable for decontamination of : <ul style="list-style-type: none"> 15L bulk bottle; buffer vessel (0.5L) along with nozzle arrangement and flexible tubing tubes, tuberlin flask, conical flask, flat round bottom flasks <p>Note:</p> <ul style="list-style-type: none"> 15L Bulk bottle- dia- 285mm, h- 430mm Tuberlin flasks (450 ml)-dia- 120 mm, h- 100 mm 250 ml conical flask-dia-85 mm, h- 140 mm 2 L Buffer conical flask- dia-165 mm, h-286 mm 3 L Flat round bottom flask- dia-185 mm, h-286 mm 5L conical Flask- dia- 220mm, h - 337mm 				
3.1.2 Articles for decontamination will be loaded manually in the autoclave so that all articles can come in contact of the steam using movable/ rotary carriage or any other better option.				
3.1.3 SS316 L loading carriage with a pair of SS316 L railing (provided inside the chamber) for smooth and easy loading. The railing should be fixed type properly welded.				
3.1.4 The chamber trolley should be provided with removable shelves for more flexibility.				
3.1.5 The chamber floor shall be on the same plane with the floor of the trolley, so that the loading carriage can directly be moved into the chamber. Loading level shall be defined by the vendor.				
3.1.6 Material will be packed in tyvek bags before loading in the equipment for decontamination.				
3.1.7 Loading environment for BF-DAT01 : loading will be from room Class D (ISO 8) cleanliness zone Loading environment for FG-DAT01 : loading will be from room Class B (ISO 5) cleanliness zone				
3.2 Brief Process Steps				
Decontamination shall have following steps <ol style="list-style-type: none"> Loading Initial Vacuum Pulsation Heating (Steaming) Hold period (Decontamination) Post vacuum Drying Unloading 				
File Name	NPI_110729_EQP_URS_BF-DAT 01_FG- DAT 01.doc		Page No.	Page 7 of 21


HLL LIFECARE LIMITED, Chennai				
Revival of BCG Vaccine Laboratory, Guindy, Chennai				
	User Requirement Specifications			
	Equipment/System	Decontamination - Autoclave		
	Identification	BF-DAT 01 & FG-DAT 01	Document	URS/BF-DAT01_FG-DAT 01
	Effective Date	2013.03.04	Revision	02
Specifications				Remarks
3.3 Output & Discharging method				
3.3.1 Articles will be unloaded from unloading side. Articles that can be reusable shall be unloaded in a class D room and taken into wash area by carriage for preparation and another round of decontamination 3.3.2 Carriage will be taken out and articles will be unloaded from the carriage. 3.3.3 All condensates and liquids shall lead to common drain without cooling. 3.3.4 All condensates and liquids from the chamber shall lead to the kill tank without cooling.				
4.0 Productivity Requirement				
4.1 Desired/ suggested capacity				
Dimension (W X H X D) and Volume (L) of autoclave chamber: BF-DAT 01: 600mm x 600mm x 600 mm, Chamber total volume: 216 L (approx) FG-DAT01: 600mm x 600mm x 600 mm, Chamber total volume: 216 L (approx)				
4.2 Standard batch size				
Not applicable				
4.3 Change Over Time				
Not applicable				
4.4 Other Productivity Requirement				
Total cycle must not exceed 2 hours.				
5.0 Containment				
Not Applicable				
6.0 GMP requirements				
6.1 Process control				
6.1.1 The equipment must operate and control the following process cycle:				
<ul style="list-style-type: none"> Vacuum leak test cycle (As per HTM 2010) 				
<ul style="list-style-type: none"> Bowie Dick cycle (17 min at 121°C and 3.5 min at 135°C) 				
<ul style="list-style-type: none"> Standard decontamination cycle (loading → steaming → hold period → slow/fast exhaust (for fluid cycle, the exhaust will be slow) 				
<ul style="list-style-type: none"> High-pressure high vacuum decontamination cycle (loading → steam/vacuum pulsing → heat up → hold period → exhaust → vacuum drying → vacuum bleeding by sterile air. 				
File Name	NPI_110729_EQP_URS_BF-DAT 01_FG- DAT 01.doc		Page No.	Page 8 of 21

HLL LIFECARE LIMITED, Chennai

Revival of BCG Vaccine Laboratory, Guindy, Chennai


HLL pharma plan®	User Requirement Specifications				
	Equipment/System	Decontamination - Autoclave			
	Identification	BF-DAT 01 & FG-DAT 01	Document	URS/BF-DAT01_FG-DAT 01	
	Effective Date	2013.03.04	Revision	02	

Specifications	Remarks
6.1.2 For the above processes following are the critical process parameters which must be controlled by the equipment	
• Pre vacuum	
• Pre pressure	
• No. of Pre pulses	
• Heat up	
• Heat up hold	
• heat up control band	
• Small valve set point	
• Decontamination hold temperature	
• Decontamination hold time	
• Temperature control band	
• Overshoot temperature	
• Decontamination stop temperature	
• Decontamination reset temperature	
• Post vacuum start pressure	
• Post vacuum	
• Post vacuum hold time	
• Post pressure	
• No of post pulses	
• Exhaust on	
• Exhaust off	
• Process end pressure	
• Chamber pressure high	
• Too long time for pre vacuum	
• Too long time for heat up	

HLL LIFECARE LIMITED, Chennai					
Revival of BCG Vaccine Laboratory, Guindy, Chennai					
HLL pharmaplan®	User Requirement Specifications				
	Equipment/System	Decontamination - Autoclave			
	Identification	BF-DAT 01 & FG-DAT 01	Document	URS/BF-DAT01_FG-DAT 01	
	Effective Date	2013.03.04	Revision	02	
Specifications					Remarks
6.2 Failure mode detection					
The Autoclave shall be capable to detect the following failure, notify the operator with alarm and shutdown the process:					
6.2.1.1 If chamber vacuum leak test is failed					
6.2.1.2 If the chamber temperature overshoots					
6.2.1.3 If chamber temperature falls below specified level & the timer stops counting					
6.2.1.4 If chamber temperature falls further below specified level & the timer resets previously counted time					
6.2.1.5 If chamber pressure is greater than the set value					
6.2.1.6 Too long time for heat up					
6.2.1.7 Too long time for pre vacuum					
6.2.1.8 Too long time for post vacuum					
6.2.1.9 If vacuum pump trips					
6.2.1.10 Door pre condition fails					
6.2.1.11 Failure in utility supply					
a) Compressed air pressure low					
b) Plant steam pressure low					
c) Softened water pressure low					
6.2.1.12 Failure in data communication					
6.2.1.13 Vendor shall propose detail list of alarms and interlocks in Functional specifications. The alarms and interlocks list shall be finalized with the final user during discussion of detail engineering design of the equipment					
6.2.1.14 Emergency stop activated					
6.2.1.15 Power failure					
6.2.1.16 Following condition need only notification to operator for procedural control					
a) UPS power low (optional)					
b) End of cycle					
c) Door opening after end of cycle					
File Name	NPI_110729_EQP_URS_BF-DAT 01_FG- DAT 01.doc			Page No.	Page 10 of 21

HLL LIFECARE LIMITED, Chennai

Revival of BCG Vaccine Laboratory, Guindy, Chennai

HLL Pharmaplan®	User Requirement Specifications				
	Equipment/System	Decontamination - Autoclave			
	Identification	BF-DAT 01 & FG-DAT 01	Document	URS/BF-DAT01_FG-DAT 01	
	Effective Date	2013.03.04	Revision	02	

Specifications	Remarks
----------------	---------


6.3 In –Process control

Manual diaphragm valves to be provided as sampling valves for steam & chamber condensate sampling.
Manual ball valve for side pocket & recirculation sampling valve
All necessary ports for steam quality testing as per EN 285 shall be incorporated.

6.4 Level of instrumentation


Sufficient and suitable instrumentation for the process, safety and productivity control as indicated in the following table:

Type of control	Purpose/ Observation	Operation range	Desired Least Count	Extent of Instrumentation			
				Indication	Alarm	Control	Recording
Temperature® (multipoint), min 5 Nos	Chamber temperatu re	0°C to + 150°C	0.1°C	Y	Y	N	Y
Temperature	Chamber condensat e drain	0°C to + 150°C	0.1°C	Y	Y	Y	Y
Temperature	Jacket temperatu re	0°C to + 150°C	0.1°C	Y	Y	Y	Y
Time	Decontam ination time	On real time basis	1 Sec	Y	Y	Y	Y
Pressure	Chamber pressure	Full vacuum to 2000 mbar	1.0 mbar	Y	Y	Y	Y
Pressure	Jacket pressure	0 to 5.0 bar	0.1 bar	Y	Y	Y	N
Differential pressure	Differentia l pressure across the sterilizing grade vacuum break filter	0 to 2000 mbar	1.0 mbar	Y	Y	N	N
Pressure	Main compress ed air line for pneumatic control	0 to 10.0 bar	0.1bar	Y	Y	N	N

HLL LIFECARE LIMITED, Chennai										
Revival of BCG Vaccine Laboratory, Guindy, Chennai										
HLL pharma plan®	User Requirement Specifications									
	Equipment/System	Decontamination - Autoclave								
	Identification	BF-DAT 01 & FG-DAT 01		Document	URS/BF-DAT01_FG-DAT 01					
	Effective Date	2013.03.04		Revision	02					
Specifications										Remarks
Pressure regulating valve along with Pressure gauge	Main plant steam line for regulating the pressure	0 to 10.0 bar	0.1bar	Y	Y	Y	N			
6.5 Batch data display and record printing										
Refer IRS (Installation Requirement Specification and Specific Instruction)										
6.6 GMP requirements (Others)										
6.6.1	Minimum 2-validation port for inserting at least 16 probes through each port during validation.									
6.6.2	Automatic F0 value calculation for each temperature monitoring location.									
6.6.3	Standard door interlocking function during decontamination cycle and at the end. <ul style="list-style-type: none">Both doors shall not be opened at a time.During the running cycle the door shall not openAfter decontamination completion the loading side door shall not be opened.After the command for unloading completion by the operator from the non sterile side, the door from loading side can be opened.									
6.6.4	Temperature trend chart recording and printing software to be provided.									
6.6.5	Vacuum pump to be provided with the system.									
6.6.6	Vacuum bleed filter: hydrophobic with arrangements for in place decontamination and provision for in-place integrity test.									
6.6.7	Provision for air leak probe as per HTM 2010									
6.6.8	Jacket to be provided with steam trap.									
6.6.9	Sampling valve in the steam inlet line for collection of steam sample.									
6.6.10	For easy & safety operation vendor shall provide the condenser in the steam sample valve outlet									
6.6.11	Sampling valve in the condensate drain line for collection of condensate sample.									
6.6.12	Vendor to give code numbers for each component.									
6.6.13	Equipment, valves, and instrumentation shall be uniquely identified in accordance with a standard numbering and location system. The system will be agreed between Vendor and Client at the time of order.									
File Name	NPI_110729_EQP_URS_BF-DAT 01_FG- DAT 01.doc				Page No.	Page 12 of 21				

HLL LIFECARE LIMITED, Chennai

Revival of BCG Vaccine Laboratory, Guindy, Chennai


nne pharmaplan®	User Requirement Specifications				
	Equipment/System	Decontamination - Autoclave			
	Identification	BF-DAT 01 & FG-DAT 01	Document	URS/BF-DAT01_FG-DAT 01	
	Effective Date	2013.03.04	Revision	02	

Specifications	Remarks
6.6.14 All valves and instruments are to be physically labeled with their equipment numbers	
6.6.15 SS panel to be flushed appropriately to the wall /ceiling/floor/LAF accordingly to avoid any dead space along with the coving on all the sides and corners	
6.7 Specific requirements	
6.7.1 Indication of chamber pressure by pressure gauge and visual LED for door open/close mounted on both non-sterile and sterile side	
6.7.2 For BF-DAT 01 Double door with vertical sliding For FG-DAT 01 Double door with vertical sliding	
6.7.3 The chamber trolleys should be provided with removable shelves for more flexibility if required. The trolley shall be 2 in numbers. The top frame is on four heavy studs for level adjustment. The rails on the top frame match with the rails in the chambers. The trolley is also provided with two fixed and two swiveling castor wheels.	
6.7.4 Fully automatic PLC/ PC based operation. Preferred make of PLC is Allan Bradley/Siemens	
6.7.5 Arrangement of alternative power supply (UPS) to control and monitoring system.	
6.7.6 The Vendor shall ensure maintenance parts availability for a minimum of 15 months from delivery.	
6.7.7 Cables, air tubes, etc required from the point (single utility point) to equipment is in scope of vendor	
6.7.8 Analogue module with back up	
6.7.9 Automatic F0 value calculation for each temperature monitoring port	
6.7.10 The chamber floor shall be on the same plane with the floor of the trolley, so that the loading carriage can directly be moved into the chamber	
6.7.11 The trolley should carry two different carriages at a time and the chamber shall also accommodate two carriages.	
6.7.12 Decontamination Chamber: The chamber shall be rectangular, with smooth and rounded corners. The chamber shall be designed as per ASME pressure vessel code. The chamber shall be made of SS316L with surface roughness less than 0.5µm. The chamber shall be re-inforced with an SS 304 jacket. The autoclave support frame for the entire structure shall be made of SS 304. The autoclave shall be able to reach and maintain decontamination temperature of 121°+ 2°C. The temperature shall be settable parameter.	

HLL LIFECARE LIMITED, Chennai				
Revival of BCG Vaccine Laboratory, Guindy, Chennai				
HLL Pharmaplan®	User Requirement Specifications			
	Equipment/System	Decontamination - Autoclave		
	Identification	BF-DAT 01 & FG-DAT 01	Document	URS/BF-DAT01_FG-DAT 01
	Effective Date	2013.03.04	Revision	02
Specifications				Remarks
6.7.13 Chamber Doors: Steam Autoclaves shall have sliding double door with automatic closing and opening. The door shall be made of SS 316L with internal surface roughness less than 0.8µm. The door gaskets shall be made of high temperature resistant silicone rubber with rounded corners				
6.7.14 Door Safety The following door safety features shall be provided for operator safety: Door interlocking to prevent simultaneous opening of both the doors. Door Process Lock to prevent opening of doors when the process is on				
6.7.15 Validation port: The chamber shall be provided with two validation ports with tri-clamp connections and with special leak tight ferrules for insertion of 16 flexible temperature sensor				
6.7.16 Vacuum Break Filter: A 0.2-micron vacuum break filter shall be provided on the sterile side for pressure equalization after vacuum creation				
6.7.17 Vendor shall install the emergency stop function on both Loading and unloading side of the decontamination autoclave				
6.7.18 Vendors have to ensure the better data restoration for both batch data recording and alarm printing. Temperature trend chart recording and printing software to be provided.				
6.7.19 Maintenance and utility shall be on one side only.				
6.7.20 Connection to drains shall be in vendor scope.				
6.7.21 If technical area is less, then there shall be a provision to remove the control panel from the machine and keep it on side wall. So that the maintenance area will increase.				
6.7.22 Equipment should be flushed with wall on both non-sterile and sterile side with bio seal.				
7.0 Constraints				
7.1 Equipment location and available space				
This equipment will be installed in the area as follows. The equipment location is indicated in the relevant block of the layout enclosed as URS Annex-1 Plant: Revival of BCG Vaccine Laboratory, Guindy, Chennai Equipment Location: BF- DAT 01				
File Name	NPI_110729_EQP_URS_BF-DAT 01_FG- DAT 01.doc		Page No.	Page 14 of 21

HLL LIFECARE LIMITED, Chennai


Revival of BCG Vaccine Laboratory, Guindy, Chennai

nne pharmaplan®	User Requirement Specifications				
	Equipment/System	Decontamination - Autoclave			
	Identification	BF-DAT 01 & FG-DAT 01	Document	URS/BF-DAT01_FG-DAT 01	
	Effective Date	2013.03.04	Revision	02	

Specifications	Remarks
<p>Floor: First Floor - Bulk</p> <p>Physical condition of the rooms:</p> <p><u>Decontamination</u></p> <p>Room dimension: 15.75 m²</p> <p>Room Height: 3000 mm</p> <p>Technical area: 4200 x 1130 mm</p> <ol style="list-style-type: none"> Room may be hazardous Class: EU Class "D" Differential Pressure: -40Pa Absolute Temperature maintained: 22°C ±2°C Relative Humidity: <55% RH <p><u>Inactivation</u></p> <p>Room dimension: 15.12 m²</p> <p>Room Height: 3000 mm</p> <ol style="list-style-type: none"> Room will be non-hazardous Class: EU Class "D" Differential Pressure: 5Pa Absolute Temperature maintained: 22°C ±2°C Relative Humidity: <55% RH <p>Equipment Location:</p> <p>FG- DAT 01</p> <p>Floor: <u>Ground Floor - Formulation</u></p> <p>Physical condition of the rooms:</p> <p><u>Air Lock</u></p> <p>Room dimension: 6.5 m²</p> <p>Room Height: 3000 mm</p> <p>Technical area: 1130 x 2420 mm</p> <ol style="list-style-type: none"> Room may be hazardous Class: EU Class "B" Differential Pressure: 35Pa Absolute Temperature maintained: 22°C ±2°C Relative Humidity: <55% RH <p><u>Washing</u></p> <p>Room dimension: 9.13 m²</p> <p>Room Height: 3000 mm</p> <ol style="list-style-type: none"> Room will be non-hazardous Class: EU Class "D" Differential Pressure: 5Pa Absolute Temperature maintained: 22°C ±2°C Relative Humidity: <55% RH 	

HLL LIFECARE LIMITED, Chennai

Revival of BCG Vaccine Laboratory, Guindy, Chennai

nne pharmaplan®	User Requirement Specifications				
	Equipment/System	Decontamination - Autoclave			
	Identification	BF-DAT 01 & FG-DAT 01	Document	URS/BF-DAT01_FG-DAT 01	
	Effective Date	2013.03.04	Revision	02	

Specifications	Remarks
----------------	---------

7.2 Utility requirement


- Electricity: 3 kW (for each Autoclave) (Report Requirement)
- Plant steam: 3-3.5 bars (Report Requirement)
- Chilled Water/Softened water: Supply: 6-7degC, Return: 11-12deg C / 3-4 bar , Ambient (Report Requirement)
- Compressed air / nitrogen pressure 6 bar g (Report Requirement)

Note: Vacuum system to be supplied by the Vendor

The vendor should plan accordingly for any change

HLL LIFECARE LIMITED, Chennai

Revival of BCG Vaccine Laboratory, Guindy, Chennai

nne pharmaplan®	User Requirement Specifications				
	Equipment/System	Decontamination - Autoclave			
	Identification	BF-DAT 01 & FG-DAT 01	Document	URS/BF-DAT01_FG-DAT 01	
	Effective Date	2013.03.04	Revision	02	


8.0 Abbreviation

List of abbreviations

HTM	Health Technical Memorandum
ISO	International Standard Organisation
LAF	Laminar Air Flow
PLC	Programmable Logic Controller
NNE	Novo Nordisk Engineering
DAT	Decontamination Autoclave
SS	Stainless steel
URS	Users requirement specification

HLL LIFECARE LIMITED, Chennai


Revival of BCG Vaccine Laboratory, Guindy, Chennai

nne pharmaplan®	User Requirement Specifications				
	Equipment/System	Decontamination - Autoclave			
	Identification	BF-DAT 01 & FG-DAT 01	Document	URS/BF-DAT01_FG-DAT 01	
	Effective Date	2013.03.04	Revision	02	

URS Annexure 1: LAYOUT POSITION


BF-DAT 01

BF026 (Decontamination) & BF046 (Inactivation Side)

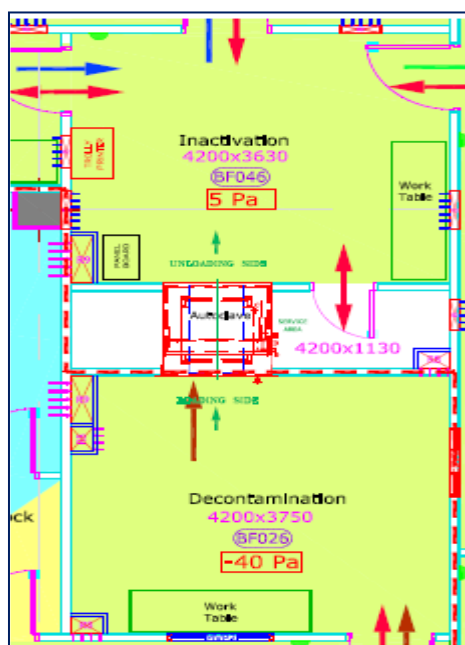

FG-DAT 01


FG037 (Airlock) & FG027 (Washing Side)

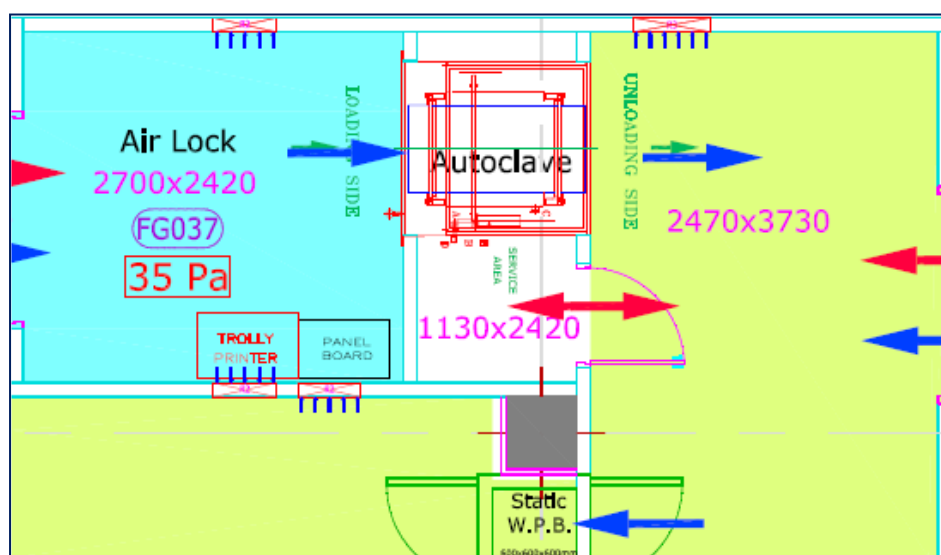

URS Annexure 2: List of Preferred Make of components

HLL LIFECARE LIMITED, Chennai


Revival of BCG Vaccine Laboratory, Guindy, Chennai

nne pharmaplan®	User Requirement Specifications				
	Equipment/System	Decontamination - Autoclave			
	Identification	BF-DAT 01 & FG-DAT 01	Document	URS/BF-DAT01_FG-DAT 01	
	Effective Date	2013.03.04	Revision	02	

SL.NO	DESCRIPTION	MAKE
A	INSTRUMENTATION	
1.	PLC	Mitsubishi/ Allen Bradley/ Siemens
2.	Operator Interface	Mitsubishi/ Allen Bradley/ Siemens
3.	Temperature transmitter	Radix/ Yokogawa/ Emerson
4.	Pressure transmitter	Siemens/ Jumo/ Wika
5.	RTD sensors	Radix/ Wika/ Waaree Instruments
6.	Temperature indicator controller	Radix/ Wika/ Waaree Instruments
7.	Printer	Epson/ HP/ Canon
8.	DC source	Shavision/ Yokogawa/ Emerson
9.	Photocell sensor	P & F/ Optex/ Metler
B	MECHANICAL	
1.	Automatic Angle Valve	Crane/ Saunder/ Gemu
2.	Manual Ball Valve	President/ Modentic/ Fluidine
3.	Needle Valve	President/ Modentic/ Fluidine
4.	Safety Valve	Teleflo/Herose/ Ciprani Harrison
5.	Non Return Valve	Leader/ Modentic/ Alfa Laval
6.	Pressure regulating valve	Klinger/ Forbes Marshall/ Armstrong International
7.	Gauges	Forbes Marshall/ Wika/ Waaree Instruments
8.	Pressure & Vacuum Switch	Orion/ Wika/ Emerson
9.	Level Switch	Mahalaxmi/ Endress & Hauser/ Emerson
10.	Steam Trap	Spirax/ Steriflow/ ITT
11.	Vacuum Break Filter	Sartorius/ Pall/ Millipore

HLL LIFECARE LIMITED, Chennai


Revival of BCG Vaccine Laboratory, Guindy, Chennai

nne pharmaplan®	User Requirement Specifications				
	Equipment/System	Decontamination - Autoclave			
	Identification	BF-DAT 01 & FG-DAT 01	Document	URS/BF-DAT01_FG-DAT 01	
	Effective Date	2013.03.04	Revision	02	

SL.NO	DESCRIPTION	MAKE
12.	Vacuum Pump	Newgenre/PPI/ Falcon Pumps
C	PNEUMATIC	
1.	Pneumatic door operating cylinder	Janatics/Rotex/ Parker
2.	Solenoid valves for door	Janatics/ Festo/ Parker
3.	Solenoid valves for Gasket	Patcon/ Festo/ Danfoss
4.	Solenoid valves for Process Valves	Janatics/ Festo/ Emerson
5.	Filter Regulator Lubricator	Janatics/ Festo/ Ingersoll
D	ELECTRICAL	
1.	Control Panel	Allen Bradley/Siemens/ ABB
2.	Limit switches	Bohmen/Siemens/ Emersen
3.	Switch gear and Relays	Siemens/ L&T/ Schneider
4.	Miniature circuit breaker	Siemens/ Havells/ Legrand
5.	Rotary switch	L&T/ Siemens/ Schneider
6.	Indication lamps	Technik / Mimic/ Schneider

HLL LIFECARE LIMITED, Chennai

Revival of BCG Vaccine Laboratory, Guindy, Chennai

nne pharmaplan®	User Requirement Specifications				
	Equipment/System	Decontamination - Autoclave			
	Identification	BF-DAT 01 & FG-DAT 01	Document	URS/BF-DAT01_FG-DAT 01	
	Effective Date	2013.03.04	Revision	02	

URS Annexure 3: List of recommended spares

SL.NO	DESCRIPTION	No. of quantity
1.	Door gasket	1 set
2.	Pressure switch	2 Nos
3.	Vacuum switch	2 Nos
4.	RTD sensors	4 No
5.	Pressure Gauge	1 Nos
6.	Compound Gauge	1 No
7.	Coil for Solenoid valves for doors	3 Nos
8.	Push buttons	2 Nos
9.	Coil for Solenoid valves for process	3 Nos
10.	Rotary Switch	2 Nos
11.	Indication bulb sets	3 Nos
12.	Temperature transmitter	2 Nos